

## MAMMOMAT 300/3000

**RX**

**Object-table top**

**Repair Instructions**

RX B7-120.091.02.01.02

Replaces:

english

06.95

	Page	Rev
<b>1 Prerequisites</b>	<b>1-1</b>	01
Tools required	1-1	
Validity	1-1	
Parts required	1-1	
<b>2 Exchange of table top</b>	<b>2-1</b>	01
Removal	2-1	
Installation	2-2	



---

## Tools required

- Standard assembly tools.

## Validity

The repair instructions are valid for:

<b>Object table</b>	<b>Part No.</b>
18x24 with grid	62 78 506
24x30 with grid	62 78 514
18x24 without grid	62 78 522
24x30 without grid	62 78 530

## Parts required

To exchange the object table top, the following parts are required:


<b>Part</b>	<b>Part No.</b>
Table top 18x24	63 14 731
or	
Table top 24x30	63 14 749



## Removal

See *Figure 1*.

- 1 Remove eight cross-section screws (1) holding the table top (2).
- 2 Remove table top (2).

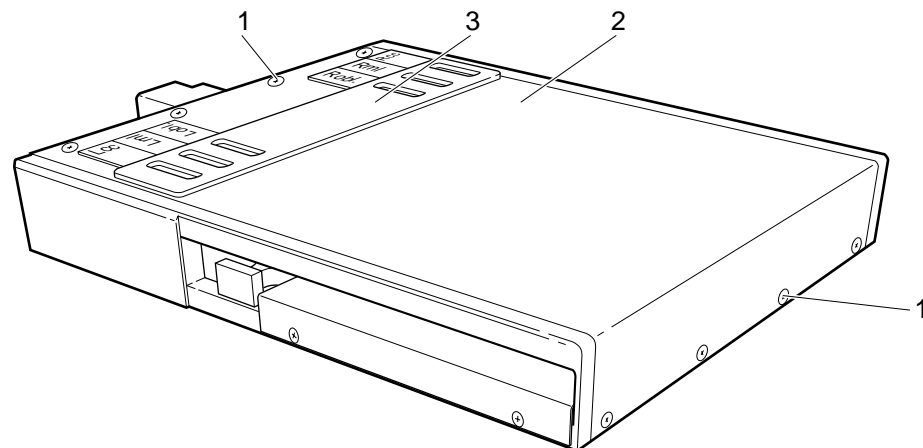


9504 0038

*Figure 1*

See *Figure 2*

- 3 Remove film marking device (3) in *Figure 1* and contact metal strip (1) by unscrewing nuts (2) and washers (3).

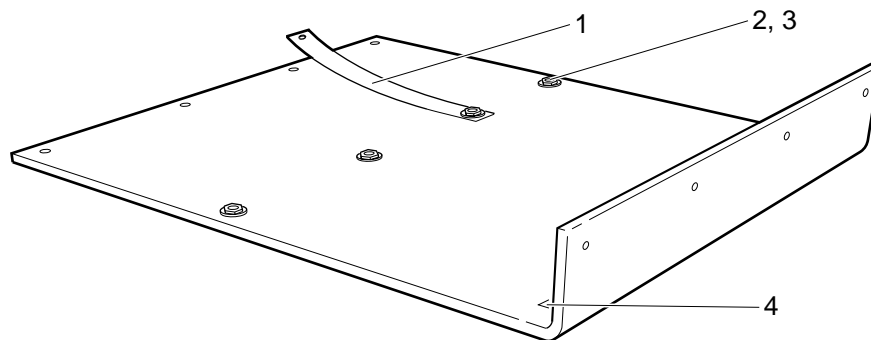

9504 0039

*Figure 2*


## Installation

- 1 Install film marking device and contact metal strip on the new table top.

**WARNING!**

*Make sure that the contact metal strip is placed in its original position (surface paint has been removed to ensure proper contact).*

- 2 Install the new table top with eight cross-section screws. The four shorter screws shall be used on the front side of the object table.
- 3 Check that the grid can be slid horizontally on the beam (4) see *Figure 2* inside the table top (if object table with grid).
- 4 Check that the film marking device is placed according to *Figure 3*.

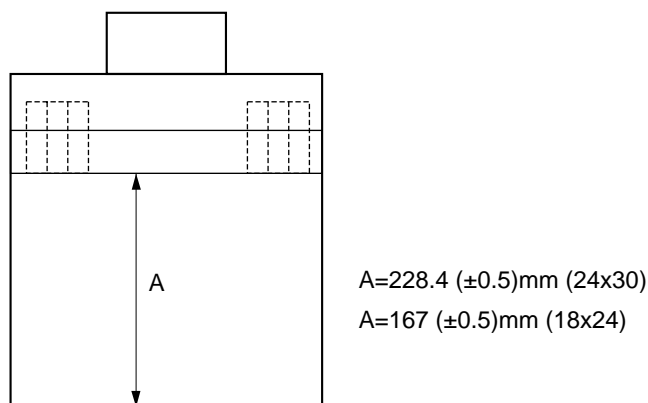

9504 0034

*Figure 3*

- 5 Fill in the LINA card and return it to the address printed on the card.